MANUXELECT



The examination and treatment couch par excellence for a small treatment room and a limited budget. Because of the scissors mechanism the ManuXelect is raised fully vertically. Consequently the ManuXelect takes up only 1.4 m2! The easily adjustable sections make this couch easy to operate and stable throughout. A cost-effective variant for someone who has little space available. The ManuXelect treatment couch comes in two 2-section models and in a 3-section variant.





- Electrical height adjustment from 48 to 105 cm
- Section adjustment with gas springs
- Ergonomic handles for operation with one hand
- Scissors movement saves space. The couch travels up and down linearly
- Padded sections with soft edges
- The nose hole (with optional plug) ensures patient comfort even when lying face down



MANUXELECT




MANUXELECT, 3-SECTION







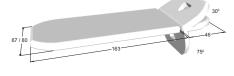

MANUXELECT, 2-SECTION, BACK SECTION








MANUXELECT




TECHNICAL SPECIFICATIONS

Size tabletop (lxw): 207 - 209 x 67 cm (optional 80 cm width)
Height adjustment: approximately 48,5-105,5 cm
Lifting time (min.-max.): approximately 35 sec.
Lifting capacity: 150 kg (optional 250 kg)
Power supply: 120/230VAC 50/60 Hz
Current consumption: 1.5 A-350W max. (light) or 2.0 A-470W max. (heavy)

Colour code frame: grey = RAL 7035 Foam density: 38 Kg/m3

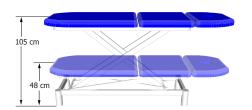

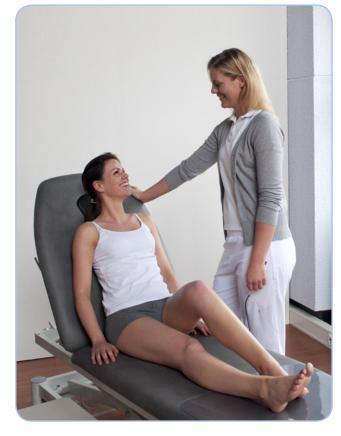

Manuxelect | | |





EXECUTIVE



012(1) black



213 grey



199(1) lavender



204 sirius blue



169(1) cornflower

blue



170(1)

baby blue

Note! Please order together with the couch because of modification of the couch






FRAME COLOURS



Ral 7035 grey

UPHOLSTERY OPTION		
Section Width	4 cm thick	
67 cm	S	
80 cm ^{(1) (2)}	0	
	S = Standard	

O = Optional

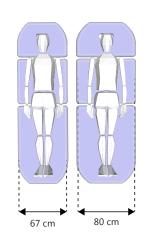

Note! Please order together with the couch because of modification of the couch (1) = Not available for traction version

(2) = Only available in colours 204, 213 and 199

A matching Enraf-Nonius couch for every interior. As in other practices, a pleasant ambience is important in physiotherapy practices. Clients must feel comfortable in the room where they are being treated. As each practice is furnished differently and has its own colour scheme, we regard it as important for the treatment couches to fit in with your look and feel. That is why we can manufacture our treatment couches in several frame and upholstery colours.

Our imitation leather upholstery contains no phthalates and is flame retardant - a unique combination that makes this material suitable for medical use!













OPTIONS MANUXELECT



HEIGHT ADJUSTMENT

- A. All-around (optional)
- B. Footswitch (optional)
- C. Handswitch (optional)







LIFTING CAPACITY

A. Electric - max. 150 kg (standard)

B. Electric - max 250 kg (optional)





OPTIONAL FEATURES

- A. Wheel raising set
- B. 1550801 Safety switch
- C. 1550804 Fold down side support (set of 2) $^{(1)}$







NOTE: Please order options together with the couch because of modification of the couch during production. $^{(1)}$ = Not suitable with 80 cm wide couches, and not suitable for traction version.





PAPER ROLL HOLDER AND PAPER ROLLS (1) 1550800 Paper roll holder for ManuXelect, for 40, 50 & 60 cm rolls, mounted under foot-section 3445280 Paper roll holder, stand-alone model 3445284 Paper roll, width 40 cm, roll-length 150 m, per 5 rolls 3445278 Paper roll, width 50 cm, roll-length 150 m, per 5 rolls 3445279 Paper roll, width 60 cm, roll-length 150 m, per 5 rolls





FIXATION ACCESSORIES (1)		
1550815	Fixation horn, black ⁽²⁾ (set of 2)	
1550806	Fixation rail for fixation horn shoulder, armpit (and for gynaecological set)	
1550807	Fixation rail for fixation horn armpit (or for gynaecological set)	
1550808	Fixation rail for fixation horn shoulder.	





FIXATION STRAPS		
3445023	Fixation strap with velcro fastener 120 cm	
3445014	Fixation strap with velcro fastener 160 cm	
3445015	Fixation strap with velcro fastener 240 cm	
3445475	Fixation strap with safety lock 250 cm	
3//5/7/	Fixation strap with safety lock 350 cm	









(2) = Requires 1550806, 1550807 or 1550808.



ACCESSORIES MANUXELECT



STOOLS (AVAILABLE IN ALL COLOURS SHOWN BELOW)3xxx010Stool, height adjustable 40-53 cm3xxx011Stool, mobile, height adjustable 59-85 cm3xxx012Stool, mobile, provided with footrest, height adjustable 59-85 cm3xxx013Stool, mobile, anatomic seat, height adjustable 43-56 cm3xxx020Stool, mobile, saddle seat, height adjustable 59-85 cm



ROLLS AND CUSHIONS (AVAILABLE IN ALL COLOURS SHOWN BELOW)				
3xxx002	Halfroll, flat bottom, small, especially for the ankle	28x9x11 cm		
3xxx003	Halfroll, flat bottom, large	56x16x18 cm		
3xxx023	Massage roll	50x Ø15 cm		
3xxx004	Neck cushion, hourglass type	30x15 cm		
3xxx005	Massage cushion	30x28x8 cm		
3xxx015	Wedge shaped cushion	50x30x15 cm		
1550813	Plug for nosehole, triangular shaped (only available in black)	16,5 x 21 cm		
3445297	Wedge shaped, Kaltenborn, rubber model (only available in black)	22,5 x 10,5 x 8,5		
3446161	Support bag with sand 1,2 kg $^{(1)}$	15x25 cm		
3446162	Support bag with sand 2 kg (1)	18x27 cm		
3446166	Support bag with sand 4 kg ⁽¹⁾	25x35 cm		
3446167	Support bag with sand 7 kg (1)	30x45 cm		





ACCESSORIES MANUXELECT



COUCH COVERS 3445240 blue Couch cover, with a cutaway for the nose 3445241 bordeaux Couch cover, with a cutaway for the nose 3445242 Couch cover, with a cutaway for the nose white 3445244 Couch cover, with a cutaway for the nose jade 3445249 Couch cover, with a cutaway for the nose dark grey 3445250 Couch cover, with a cutaway for the nose royal blue 3445251 Couch cover, with a cutaway for the nose camel 3445261 Couch cover, without a cutaway for the nose blue 3445034 Couch cover, without a cutaway for the nose bordeaux 3445267 Couch cover, without a cutaway for the nose white 3445262 Couch cover, without a cutaway for the nose 3445248 Couch cover, without a cutaway for the nose dark grey 3445247 royal blue Couch cover, without a cutaway for the nose 3445252 Couch cover, without a cutaway for the nose camel Elastic straps for couch cover with clip fittings 3445259









GYNAECOLOGY ACCESSORIES (1)			
1550806	Fixation rail for gynaecological set (and for fixation horn shoulder, armpit)		
1550807	Fixation rail for gynaecological set (and for fixation horn armpit)		
1550810	Gynaecological Comfort set ⁽²⁾		
1550805	Gynaecological set ⁽²⁾		

(1) = Not suitable for traction version or 80 cm wide couches.

(2) = Requires 1550806 or 1550807









MISCELLANEOUS

3446261 Set of 3 exercise stools, 50 x 38 x 32 cm, $40 \times 32 \times 21$ cm and $30 \times 22 \times 14$ cm













The traction version of the ManuXelect has a rolling section with a range of motion of 15 cm to perform traction therapy. This rolling section can also be locked in place for conventional use.

Your ManuXelect traction couch comes standard with a fixation rail on the head section and a traction frame for mounting the Eltrac 471 traction device.

The options for the height adjustment and the lifting capacity of the ManuXelect traction couch are the same as for the other ManuXelect versions.

Attention! The ManuXelect Traction is available only in the colours grey (colour number 213) and Sirius blue (colour number 204) and only with a section width of 67 cm.

TRACTION ACCESSORIES Flexistool grey, with separate height adjustment for 3xxx014 each leg, available in different colours 3445022 Pelvic belt for Manumed Traction 3445021 Thoracic belt for Manumed Traction Saunders clinical traction device including mounting 1471808 bracket for the Eltrac 471 1671901 Eltrac 471, unit for traction therapy











TECHNICAL SPECIFICATIONS

Size tabletop (lxw): approximately 240 x 67 cm (incl. traction frame)

