

comfort.



**MULTILINE NEXT BT AND BT PLUS. OVERCHAIR TABLES.**



**GREINER**

**BT and BT plus overchair tables:**


- ◆ Practical, sturdy and convenient overchair tables, easy to move around
- ◆ Sturdy metal frame is safe for loads up to 30 kg
- ◆ The height can be easily adjusted with one hand – by you or your patient
- ◆ Mechanism for tilting the tray by 20° or 40° on both sides (BT plus)
- ◆ Removable trays with recessed grips, securely fixed
- ◆ Very hygienic thanks to easy-clean surfaces
- ◆ Space-saving storage – the tables can be nested close together



**MULTILINE NEXT BT AND BT PLUS. OVERCHAIR TABLES FOR HOSPITAL USE.**

They roll in softly: the new GREINER BT and BT plus overchair tables. Developed and designed as top-quality additions to the MULTILINE NEXT DC and IT treatment chairs, these side tables ensure the focus is firmly on patient comfort. The BT with fixed tray and the BT plus with swivelling tray are flexible overchair tables that allow you to ensure your patients are as comfortable as possible when undergoing dialysis, chemotherapy or other treatments. They provide a unique combination of stability and lightness and can be adjusted with one hand. Everything you need to care for your patient is right at your fingertips.

**BT plus overchair table**  
Variable height adjustment with one hand, swivelling tray, hook for accessories and holder for magazines or patient files.

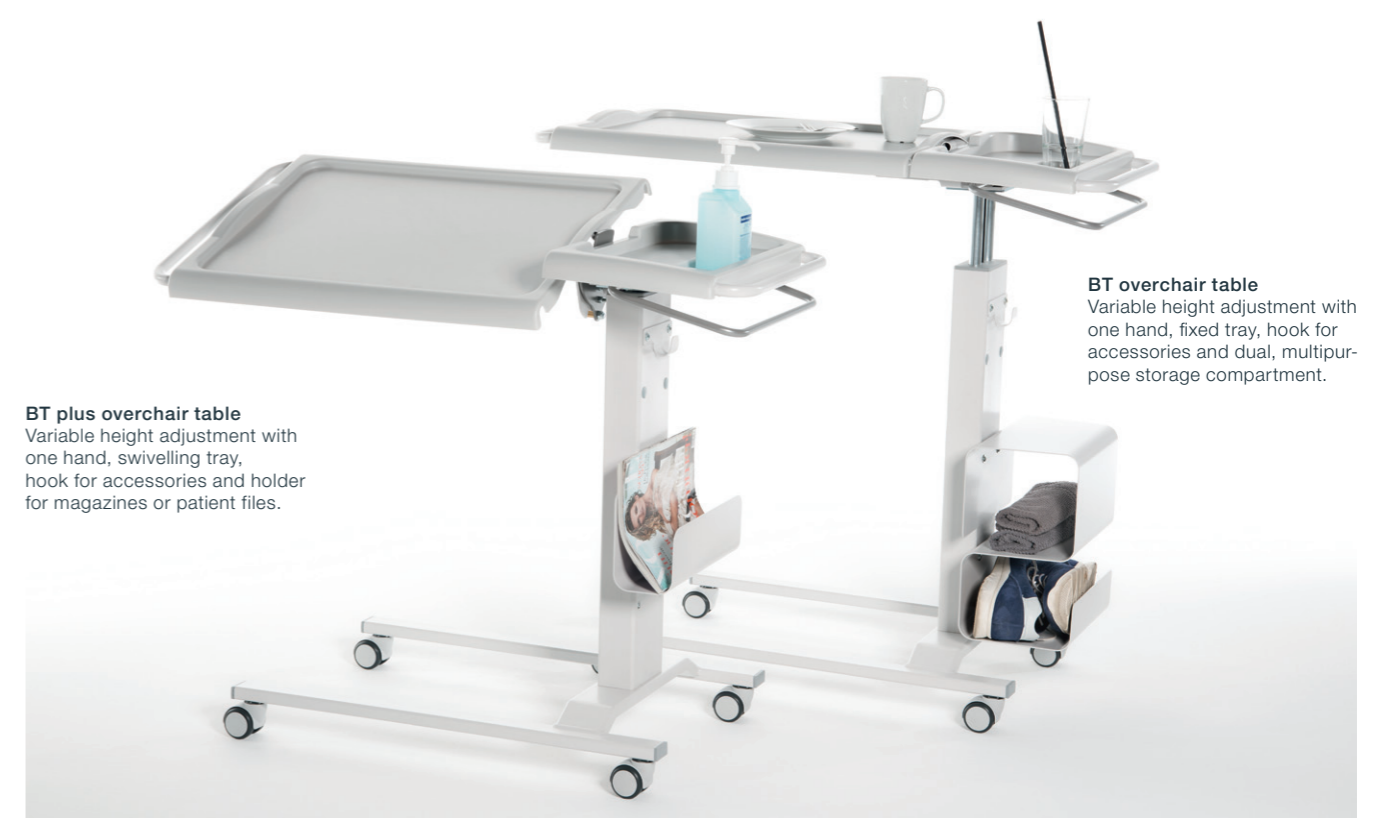

**BT overchair table**  
Variable height adjustment with one hand, fixed tray, hook for accessories and dual, multipurpose storage compartment.



Perfect for combining with the MULTILINE NEXT DC and IT (IT shown here)




Height adjustment with dual control



Swivelling tray (BT plus)



Design and quality made by GREINER

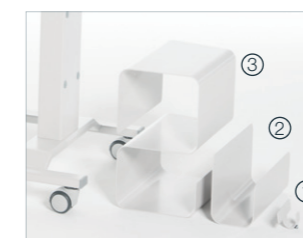

Multipurpose accessories, see table

The epitome of functionality, ergonomics, stability and hygiene, the new overchair tables have been meticulously designed to ensure they meet the most stringent requirements of patients and nurses alike. The sturdy frame in typical GREINER quality is combined with smooth-running castors to ensure the table is easy to move while still safely holding loads of up to 30 kilos. Powerful pneumatic springs make it easy to gently adjust the height of the tray with one hand – and provides added safety if the patient is in the emergency Trendelenburg position. On the BT plus the tray can be tilted at the touch of a button so that it is ergonomically positioned for both patient and nurse. The removable trays have recessed grips and are the perfect size for holding meal trays, patient items, clinical accessories and disinfectants. The tough, resistant surfaces are quick and easy to clean and disinfect. A range of accessories are also available to improve the patient experience still further.

**Dimensions, technical specifications and accessories**

|                                |                        |
|--------------------------------|------------------------|
| ◆ Overchair tables (L x W)     | 97 cm x 46 cm          |
| ◆ Tray surface (L x W)         | 56 cm x 39 cm          |
| ◆ Side tray (L x W)            | 37 cm x 21 cm          |
| ◆ Variable height adjustment   | 25 cm (79 – 104 cm)    |
| ◆ Diameter of castors          | 50 mm                  |
| ◆ Tilt angle of tray (BT plus) | 20° / 40° on each side |
| ◆ Maximum load                 | 30 kg                  |
| ◆ Weight                       | Approx. 22.5 kg        |
| ◆ Colour                       | RAL 7047               |

**Accessories**

|  |                                    |
|--|------------------------------------|
| ◆ Hook (1)                                   | For bags, etc.                     |
| ◆ Magazine holder (2)                        | For magazines, patient files, etc. |
| ◆ Dual, multipurpose storage compartment (3) | For shoes, wipes, books, etc.      |



## we design quality.

Agreeably health-promoting seating has motivated us since 1922. For decades we have been making high-quality, technically sophisticated equipment within the medbest range. Our innovative chairs for examination, dialysis, oncology, blood sampling, care and transport are all made in Germany. Our quality commitment is verified by an independently audited quality assurance system: GREINER is certified under DIN EN ISO 9001 and 13485.



Test mark for  
quality management system



GREINER GmbH  
Wettstrasse 1 · 74385 Pleidelsheim · Germany  
T +49 7144 81120 · F +49 7144 8112 99  
info@greiner-gmbh.de · www.greiner-gmbh.de

