

MANUMED OPTIMAL 2-SECTION



With the Manumed Optimal 2-section treatment couch you have the modern Manumed Optimal quality couch at your disposal. Robust in appearance and simple to operate. Stable and durable.

Less is more with the Manumed Optimal 2-section couch. The 2-part head section, the 2-part head section with separate arm supports and the 2-part back section offer a full range of options to the therapist looking for an uncomplicated couch for safe treatment.



*"Easy change from sitting to lying position.
Ideal for electrotherapy treatment."*




- 2 (up to 4) adjustable sections
- Load capacity up to 200 kg (optional 250 kg)
- Head section folds down completely for better access of the cervical region
- Optimally positioned handles
- Long nose hole, allowing proximal or distal movement of the patient even when lying in prone position
- Optional a triangular nose hole or without nose hole
- Wheel raising mechanism for easy movement of the couch
- Height can be adjusted hydraulically (no mains lead needed) or electrically (very convenient, smooth and quick)
- Completely relaxed starting position in the face-down position with separate head section with adjustable arm rests (optional)



ENRAF-NONIUS B.V.
Vareseweg 127
3047 AT Rotterdam
The Netherlands
T: +31-(0)10 - 203 06 00
E: info@enraf-nonius.nl
www.enraf-nonius.com

MANUMED OPTIMAL 2-SECTION


MANUMED OPTIMAL 2-SECTION, HEAD-SECTION

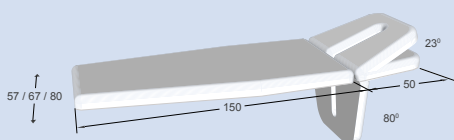

Type

- 311 2-section, head-section hydraulic H/L
- 321 2-section, head-section electric H/L


MANUMED OPTIMAL 2-SECTION, HEAD-SECTION WITH ADJUSTABLE ARM SECTIONS

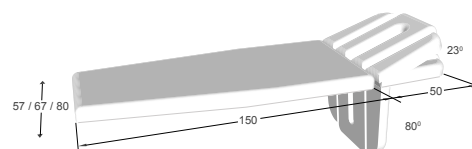

Type

- 313 2-section, head-section with adjustable arm sections, hydraulic H/L
- 323 2-section, head-section with adjustable arm sections, electric H/L

MANUMED OPTIMAL 2-SECTION, BACK SECTION

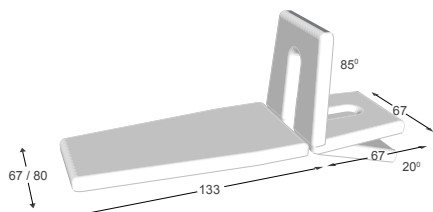

Type

- 312 2-section, back-section, hydraulic H/L
- 322 2-section, back-section, electric H/L

TECHNICAL SPECIFICATIONS

- Size tabletop (lxw): 201 x 67 cm (optional 57 or 80 cm width)
- Height adjustment: 45-95 cm (backsection)
47-97 cm (headsection)
- Lifting time (min.-max.): approximately 18 sec.
- Lifting capacity: 200 kg
(electrical model optional 250 kg)
- Power supply: 120/230VAC 50/60 Hz
- Current consumption: 2.0 A max.

NOTE: ALL-AROUND OPERATION IS AN OPTION.



MANUMED COLOURS

A matching Enraf-Nonius couch for every interior. As in other practices, a pleasant ambience is important in physiotherapy practices. Clients must feel comfortable in the room where they are being treated. As each practice is furnished differently and has its own colour scheme, we regard it as important for the treatment couches to fit in with your look and feel. That is why we can manufacture our treatment couches in several frame and upholstery colours. Our imitation leather upholstery contains no phthalates and is flame retardant – a unique combination that makes this material suitable for medical use!




NATURAL



811



850

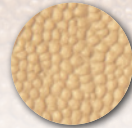


851



886



888

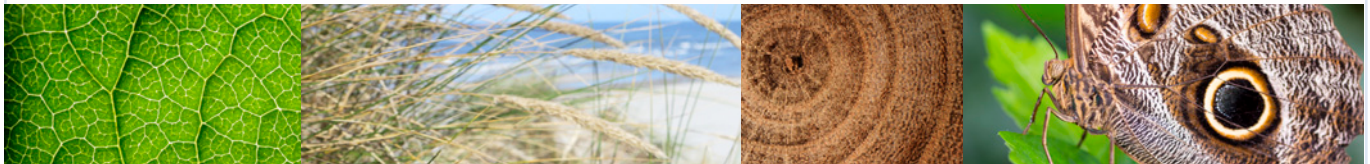

LIFESTYLE



410
anthracit



406
petrol



407
cherry



408
terracotta




409
skin




EXECUTIVE



012
black

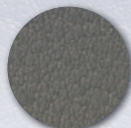


213
grey

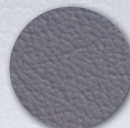


199
lavender



204
sirius blue

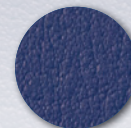

169
cornflower
blue

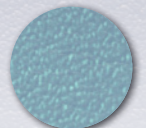

170
baby blue



Frame colours



Ral 7035
grey



Ral 9005
black




Ral 9010
white

OPTIONS⁽¹⁾ MANUMED OPTIMAL 2-SECTION

UPHOLSTERY OPTIONS

Section Width	Standard 4 cm thick	DeLuxe 6 cm thick	Extra hard 4 cm thick	Extra soft 4 cm thick	Trapzium shaped leg section	Rectangular shaped leg section
57cm	s	o	o	o	x	s
67 cm	s	o	o	o	s	o
80 cm	s	o	o	o	s	x

s = standard o = optional x = not available

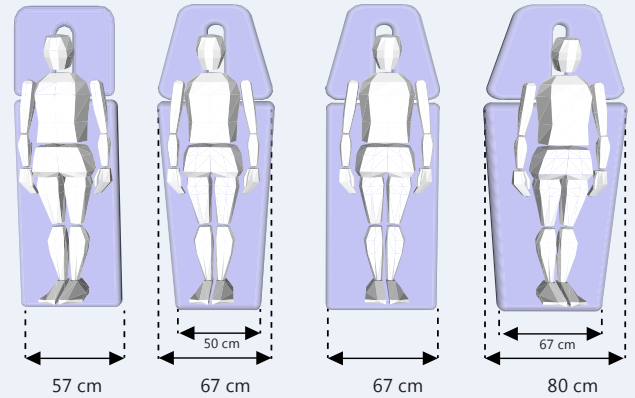


FOAM DENSITY:

Standard: 38 Kg/m³
 DeLuxe: 38 Kg/m³
 Extra hard: 70 Kg/m³
 Extra soft: 4 cm 38 Kg/m³ + 1,5 cm 30 Kg/m³



NOSE-HOLE

- A. Oblong nose-hole including plug (standard)
- B. Triangular shaped nose-hole (optional) excluding plug, please order seperately 1550813
- C. No nose-hole (optional)

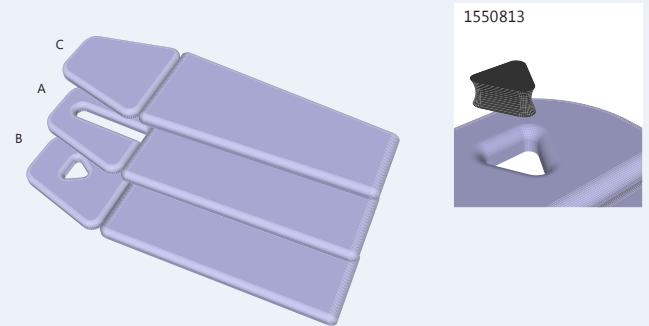

HEIGHT ADJUSTMENT

- A. Hydraulic (311 / 312 / 313)
- B. Electric (321 / 322 / 323) Footswitch (standard)
- C. Electric (321 / 322 / 323) All-around (optional)



LIFTING CAPACITY

- A. Electric - max. 200 kg (standard 321 / 322 / 323)
- B. Electric - max 250 kg (optional 321 / 322 / 323)
- C. Hydraulic - max 200 kg (311 / 312 / 313)



OPTIONAL FEATURES

- A. Safety switch)



(1) = NOTE: Please order options together with the couch because of modification of the couch during production

ACCESSORIES MANUMED OPTIMAL 2-SECTION



PAPER ROLL HOLDER AND PAPER ROLLS

3450104	Paper roll holder for Manumed Optimal, for 40, 50 & 60 cm rolls, mounted under foot-section (not for type 312, 322)
3450105	Paper roll holder for Manumed Optimal, for 40, 50 & 60 cm rolls, mounted under headsection
3450065	Paper roll holder foot-section for type 312, 322
3445280	Paper roll stand (floor type) for 40, 50 & 60 cm rolls
3445284	Paper roll, width 40 cm, roll-length 150 m, per 5 rolls
3445278	Paper roll, width 50 cm, roll-length 150 m, per 5 rolls
3445279	Paper roll, width 60 cm, roll-length 150 m, per 5 rolls



FIXATION ACCESSORIES

3xxx035	Shoulder/neck support, set of 2 (available in different colours, see previous pages)
3450036	Clamp for fixation horn 3xxx035 (2 pieces) ¹
3XXX028	Armpit fixation horn (available in different colours, see previous pages)

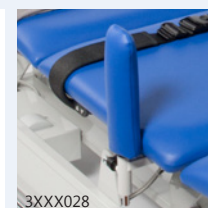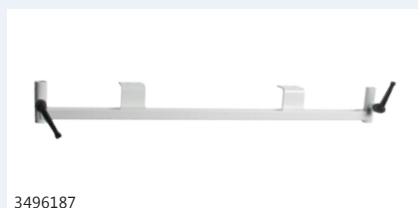

FIXATION STRAPS

3445023	Fixation strap with velcro fastener 120 cm
3445014	Fixation strap with velcro fastener 160 cm
3445015	Fixation strap with velcro fastener 240 cm
3445475	Fixation strap with safety lock 250 cm
3445474	Fixation strap with safety lock 350 cm



ACCESSORIES MANUMED OPTIMAL 2-SECTION

STOOLS (AVAILABLE IN ALL COLOURS SHOWN BELOW)

3xxx010	Stool, height adjustable 40-53 cm
3xxx011	Stool, mobile, height adjustable 59-85 cm
3xxx012	Stool, mobile, provided with footrest, height adjustable 59-85 cm
3xxx013	Stool, mobile, anatomic seat, height adjustable 43-56 cm
3xxx020	Stool, mobile, saddle seat, height adjustable 59-85 cm



ROLLS AND CUSHIONS (AVAILABLE IN ALL COLOURS SHOWN BELOW)

3xxx002	Halfroll, flat bottom, small, especially for the ankle	28x9x11 cm
3xxx003	Halfroll, flat bottom, large	56x16x18 cm
3xxx023	Massage roll	50x Ø15 cm
3xxx004	Neck cushion, hourglass type	30x15 cm
3xxx005	Massage cushion	30x28x8 cm
3xxx015	Wedge shaped cushion	50x30x15 cm
1550813	Plug for nosehole, triangular shaped (only available in black)	16,5 x 21 cm
3445297	Wedge shaped, Kaltenborn, rubber model (only available in black)	22,5 x 10,5 x 8,5
3446161	Support bag with sand 1,2 kg ⁽¹⁾	15x25 cm
3446162	Support bag with sand 2 kg ⁽¹⁾	18x27 cm
3446166	Support bag with sand 4 kg ⁽¹⁾	25x35 cm
3446167	Support bag with sand 7 kg ⁽¹⁾	30x45 cm


(1) = only available in grey




EXECUTIVE



NATURAL

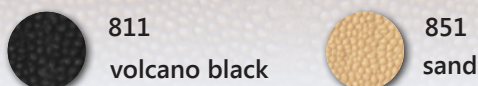


LIFESTYLE

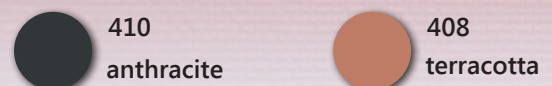

ACCESSORIES MANUMED OPTIMAL 2-SECTION

COUCH COVERS

3445240	Couch cover, with a cutaway for the nose	blue
3445241	Couch cover, with a cutaway for the nose	bordeaux
3445242	Couch cover, with a cutaway for the nose	white
3445244	Couch cover, with a cutaway for the nose	jade
3445249	Couch cover, with a cutaway for the nose	dark grey
3445250	Couch cover, with a cutaway for the nose	royal blue
3445251	Couch cover, with a cutaway for the nose	camel
3445261	Couch cover, without a cutaway for the nose	blue
3445034	Couch cover, without a cutaway for the nose	bordeaux
3445267	Couch cover, without a cutaway for the nose	white
3445262	Couch cover, without a cutaway for the nose	jade
3445248	Couch cover, without a cutaway for the nose	dark grey
3445247	Couch cover, without a cutaway for the nose	royal blue
3445252	Couch cover, without a cutaway for the nose	camel
3445259	Elastic straps for couch cover with clip fittings at both ends, set of 3	

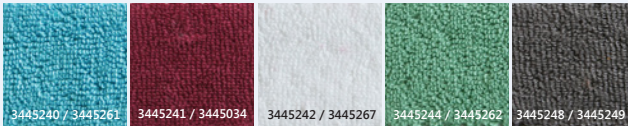

MISCELLANEOUS

3450067 Raising set (30 mm)



3450067

