

# Sonopuls 190

Ultrasound therapy:  
fast and easy

## Sonopuls 190

The new Sonopuls 190 is fast and easy. A device with an extremely fast start-up, a fast way to set-up treatment parameters and a device that gives fast therapeutic results. The compact and ergonomic design, the clear TFT screen and the touchscreen interface are making this small device a big asset for the therapist.

Ultrasound therapy is a proven effective treatment method for numerous pathologies. In order to simplify the therapy selection, a number of preprogrammed protocols are stored in the Sonopuls 190. These evidence based protocols provide you with clear instructions and will set the device parameters for treatment automatically.



- Evidence based protocols
- Programmable
- Extremely quick start-up
- Crystal clear tft colour screen
- Touch screen user interface
- Programmed for speed
- Ergonomic design
- Contact control
- 2 Treatment heads possible
- Energy efficient
- Compact and stylish
- Ultra-flexible cable
- By Enraf-Nonius: safe, reliable and durable

### Evidence based protocols

Provided with 25 evidence based treatment protocols, including anatomical representations and clear guidelines for treatment.

### Programmable

Besides the 25 evidence based protocols you can also store 20 protocols according to your own guidelines.

### Extremely quick start-up

Start the unit and it is ready for immediate use. No start-up delay, no operating system to load but instant availability of the full functionality of the unit.

### Crystal clear TFT colour screen

A crystal clear and scratch resistant screen that allows full visibility on all parameters from all viewing angles.

### Touch screen user interface


Easy to use by means of a clear menu via a touch screen interface. Extremely intuitive and as easy as a state-of-the-art smart phone.

### Programmed for speed

The menu structure is optimized to speed. In most cases you can start a pathology-specific treatment in just 3 steps! Easy, quick and simple.

### Ergonomic design

The treatment head is made to fit your hand comfortably. A design completely based on functionality and user friendliness.

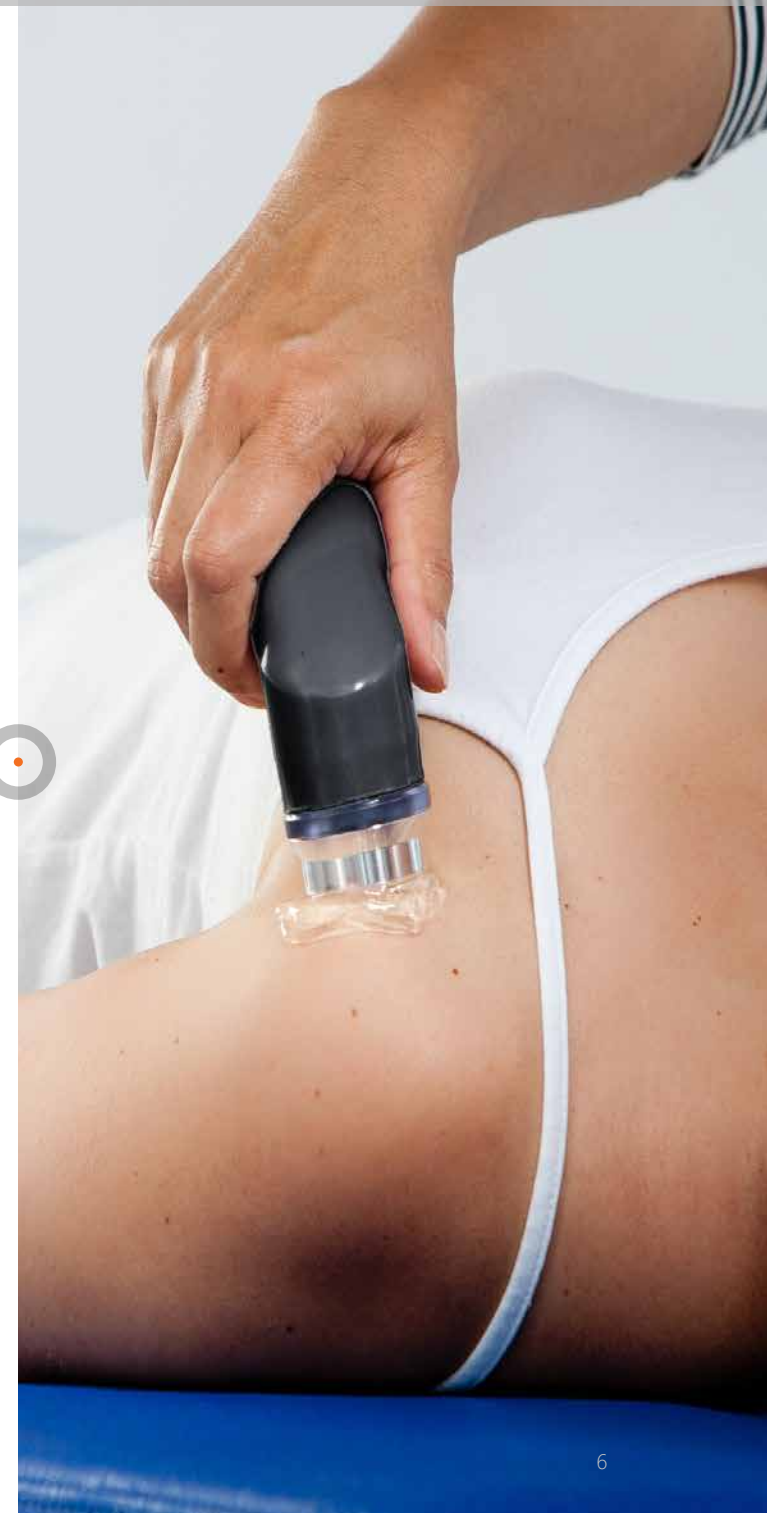



### Contact control

A 360 degrees light ring for contact control. A clear indication of the effective transfer of energy that enables you to conduct the therapy as efficient as possible.

### 2 Treatment heads possible

You can connect both the small 0,8 cm<sup>2</sup> as well as the large 5 cm<sup>2</sup> ultrasound treatment head in order to quickly change heads.

### Energy efficient

Sounds as a paradox for a unit which primary goal is to transfer energy to the human body. Due to smart energy saving modes no waste is generated when the device is in rest-mode.

### Compact and stylish

An extremely compact device, small and light-weight for easy transportation. At the same time stylish and beautiful, a device that is worth-while showing!

### By Enraf-Nonius: safe. reliable and durable

An Enraf-Nonius device, the worldwide supplier for physiotherapy and rehabilitation equipment. Made with special attention to quality, durability and safety.

### A treatment head with added value

An intelligent treatment head, packed with electronics and covered with a special coating making it smooth and comfortable to use and easy to clean.

### Ultra-flexible cable

The cable between the unit and the ultrasound head is extremely flexible. This assures ease of use. The cable adapts according to your movements and will never impede any motion of the treatment head.

# Sonopuls 190

## 1631901 Sonopuls 190

With multi-frequency ultrasound head large of 5 cm<sup>2</sup>

### Standard accessories

- 1630905 Ultrasound multifrequency treatment head large (5 cm<sup>2</sup>)
- 1631820 Head holder for multi-frequency treatment head
- 3442929 Contact-gel, 1 bottle of 250 ml <sup>(1)</sup>
- 1631750 User manual Sonopuls 190 (on CD)
- 0167154 Information sheet ultrasound gel
- 0167278 Information sheet mounting US head holder(s) Sonopuls 190
- 3440001 Screwdriver
- 3444357 Mains connecting cable



### Optional accessory

- 1630915 Ultrasound multifrequency treatment head small (0,8 cm<sup>2</sup>)



1630915

(1) = The Sonopuls is delivered with 1 bottle of contact-gel. The articlenumber 3442929 however represents a box of 12 bottles.



# Sonopuls 190

## 1631904 Sonopuls 190 with StatUS™ Pack 100

With StatUS™ applicator and applicator holder



### Standard accessories for 1631904

- 0169842 Head holder StatUS™ applicator
- 3442941 Gelpad StatUS™ (240x)
- 3442942 Fixation ring gelpad StatUS™ (3x)
- 1631750 User manual Sonopuls 190
- 1629750 User manual StatUS™ Pack 100
- 1629751 Information booklet StatUS™ Pack 100
- 3440001 Screwdriver
- 3444357 Mains connecting cable

# ORDERING DATA

## 1629901 StatUS™ Pack 100

With StatUS™ applicator and applicator holder



### Standard accessories for 1629901

- 0169842 Head holder StatUS™ applicator
- 3442941 Gelpad StatUS™ (240x)
- 3442942 Fixation ring gelpad StatUS™ (3x)
- 1629750 User manual StatUS™ Pack 100
- 1629751 Information booklet StatUS™ Pack 100



# Sonopuls 190

## TECHNICAL SPECIFICATIONS

Ultrasoundfrequencies:	1 en 3 MHz
Ultrasound:	continuous and pulsed
Duty cycles	5, 10, 20, 33, 50, 80%
Pulse rate:	16 Hz, 48 Hz and 100 Hz
Numbers of US connections:	2
Intensity	0-2 W/cm <sup>2</sup> continuous, 0-3 W/cm <sup>2</sup> pulsed <sup>(1)</sup>
ERA US treatment head large:	5 cm <sup>2</sup>
ERA US treatment head small:	0,8 cm <sup>2</sup>
ERA StatUS™ applicator:	5 cm <sup>2</sup>
Pre-programmed protocols:	25 evidence based protocols for Ultrasound therapy 52 evidence based protocols for StatUS™ therapy
Free program positions:	20 pcs
TFT touchscreen:	480x272 pixels
Main power:	100-240 V ± 10% (50/60Hz)
Dimensions (wxdxh):	21x19x9 cm
Weight:	1152 gr
This is a medical device.	It complies with all applicable medical regulations. For more details please consult the instructions for use (downloadable via <a href="http://www.enraf-nonius.com">www.enraf-nonius.com</a> ).

(1) = At a duty cycle of 80% the maximum output is 2,5 W/cm<sup>2</sup>.



# Sonopuls 190

# ACCESSORIES

## Ultrasound contact-gel

3442929 Contact-gel, bottle 250 ml, box of 12

3442931 Contact-gel, canister of 5 l

## Cases and bags

3444061 Bag for Endomed and Sonopuls 1-serie, black



3442929



3442931



3444061





## Sonopuls 190

Ultrasound therapy:  
fast and easy

Enraf-Nonius B.V. | Vareseweg 127 | 3047 AT Rotterdam | The Netherlands  
[www.enraf-nonius.com](http://www.enraf-nonius.com) | [info@enraf-nonius.nl](mailto:info@enraf-nonius.nl) | + 31 - ( 0 ) 10 20 30 600



PARTNER FOR LIFE