



## EN-Tree Pulley

A complete series of  
multi-functional pulley systems



PARTNER FOR LIFE



## EN-Tree Pulley

The Enraf-Nonius EN-Tree assortment offers a complete gamma for medical training therapy (MTT). Medical training therapy is a concept that is used for:

- strengthening muscles
- mobilisation
- improving one's endurance
- improving one's co-ordination

MTT makes use of various aids and exercise systems. In addition to our EN-Dynamic exercise equipment and EN-Cardio training equipment, we offer a broad range of pulley systems. EN-Tree pulleys are ideal for MTT because every movement can be simulated and trained. In addition to isolated movements, pulleys can be used to mimic almost every movement that one makes in daily life. The EN-Tree systems are well accessible, they enable one to gradually increase the training programme and they require very little space.

EN-Tree renders medical training accessible for everyone!

## EN-Tree Pulley

The EN-Tree products pallet comprises various devices.

First of all, we offer 4 different pulleys:

- 1 EN-Tree Pulley (24 kg), the standard versatile pulley with a range of 0.25 kg to 24 kg.
- 2 EN-Tree Pulley Explosive (54 kg), the ideal pulley for training rapid explosive movements with a range of 0.33 kg to 18 kg
- 3 EN-Tree P, the pneumatic pulley for both power training as well as for training rapid explosive movements with a range of 0.66 kg (6.6N) to 24 kg! Moreover, the EN-Tree P can be adjusted fully automatically using the EN-Train system.
- 4 EN-Tree M, the measuring pulley, available in a normal or explosive version and equipped with software for the purpose of quantifying the exercises with the pulley.

In addition, we offer a wide range of accessories for all of our pulleys. These accessories enable one to carry out very specific exercises.

We also offer an exercise bench in addition to the various pulleys. This bench is ideal for exercising with the EN-Tree pulleys. In addition, this bench is also ideal for mobilisation exercises and muscle-strengthening exercises without additional equipment. Thanks to its small size and the small wheels the EN-Tree Bench can easily be moved elsewhere.





## EN-Tree Pulley

- Multi-functional
- Accurate dosage of the exercise resistance
- Complete possibility to exercise on a surface area of 2.5m<sup>2</sup>
- Extremely smooth dual-bearing rollers
- Height adjustable with indicator
- Adjustable cord extension
- Very thin but extremely solid cord (less resistance from the rollers)


EN-Tree  
renders medical training  
accessible for everyone!

## EN-Tree P

The EN-Tree P (for pneumatic) is a pulley system that does not make use of the traditional weights, but rather a pneumatic resistance. In addition to the advantages of a traditional pulley system, such as versatility, requires little space and good accessibility, the pneumatic system like the EN-Tree P has a few more advantages:

- an extremely low minimum load
- a high maximum resistance (720 N)
- it is possible to increase the load in very small steps
- Silent during exercising
- no mass inertia due to the pneumatic resistance, making the EN-Tree P the ultimate Explosive pulley
- a pulley that can be fully controlled (and adjusted automatically) from EN-Train, enabling one to undergo consistent and controlled training sessions

With that, EN-Tree P is a unique system that offers added value to any all-round practice, with or without EN-Train.

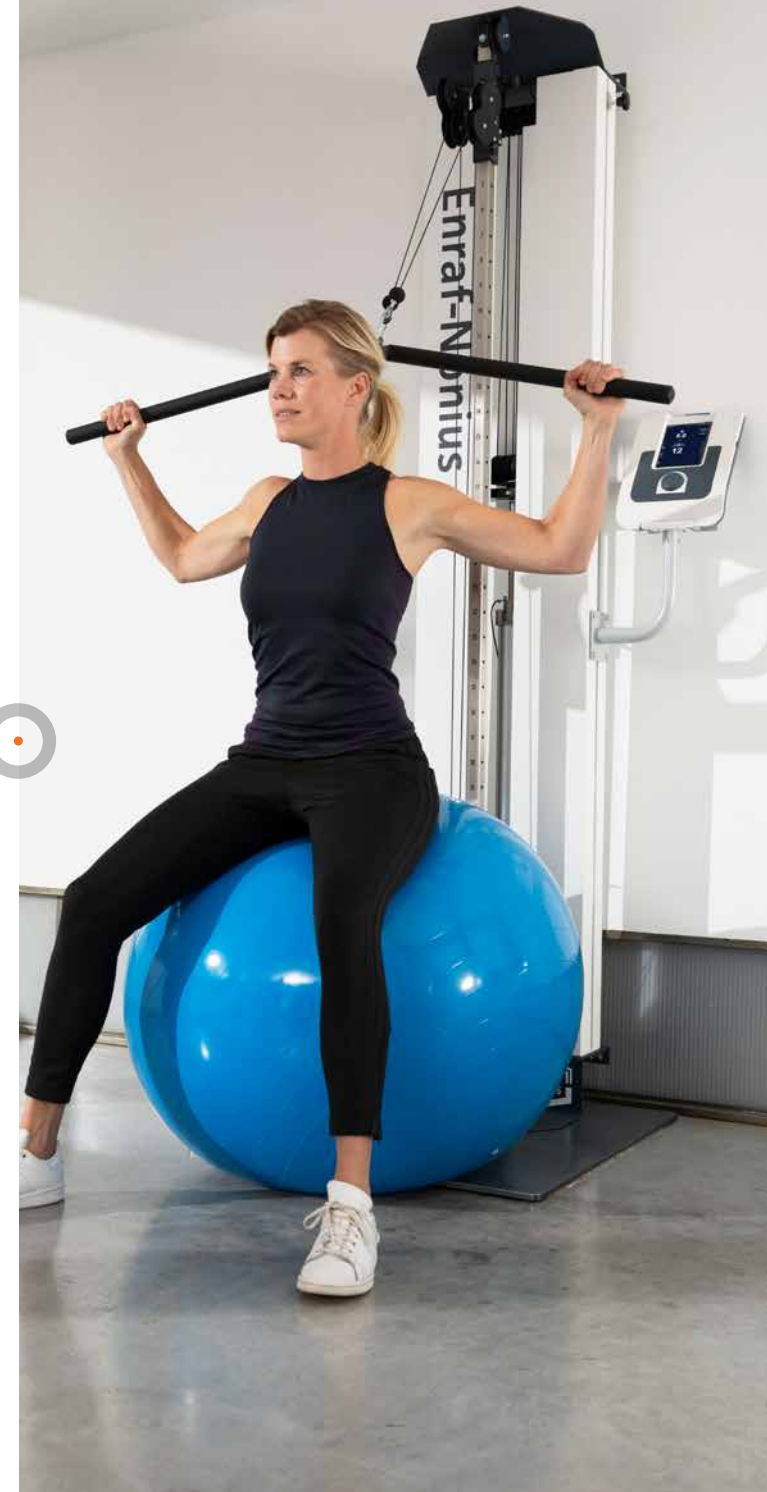

## EN-Tree M

In physical therapy and rehabilitation it is very important to quantify the functional capacity of the human muscular-skeletal system. This is in order to determine the appropriate individual load during exercise therapy or to evaluate the progress during rehabilitation. Custom systems for the measurement of functional capacity require many square meters of floor space and high investments. To make measurement techniques available for every therapist Enraf-Nonius offers the EN-Tree M, the unique unit for extremely accurate physical therapy and rehabilitation.

The EN-Tree M is the measurement edition of the well-known EN-Tree Pulley, the ideal aid for functional training and rehabilitation. The EN-Tree M accurately determines position, velocity, force, power and work. With these parameters you can quantify the functional capacity of your patients with a PC.

In the powerful software you can compose various test and training protocols. To analyse the test and training results a professional and complete package is provided which allows you to assess optimum individual treatment.





## EN-Tree M

The goals of treatment can be determined by:

- left-right comparison
- assessment of normal or abnormal movement patterns from the graphs;
- comparison of test results between individuals

Progress during rehabilitation can be monitored by comparing the results of different tests recorded in time.

For a reliable reproduction of a training session or a test you initiate the pulley set-up as well as the positioning of the patient.



## EN-Tree Pulley (24 kg)

With EN-Tree products, active exercise is possible for almost anyone. Healthy individuals, patients with specific exercise requirements and wheelchair users can enter into their own training programs. This is possible because of the wide variety of load application and the ease of adaptation of the EN-Tree to every individual's need. The EN-Tree Pulley allows an increase of the weight from 0,25 to 24 kg in small steps. The EN-Tree Pulley, allows the exercise of almost every human movement. This includes mono-articular movements as well as complex displacements.


## ORDERING DATA



1455970 EN-Tree Pulley (24 kg)




1455971 EN-Tree Pulley MD (24 kg with cover)

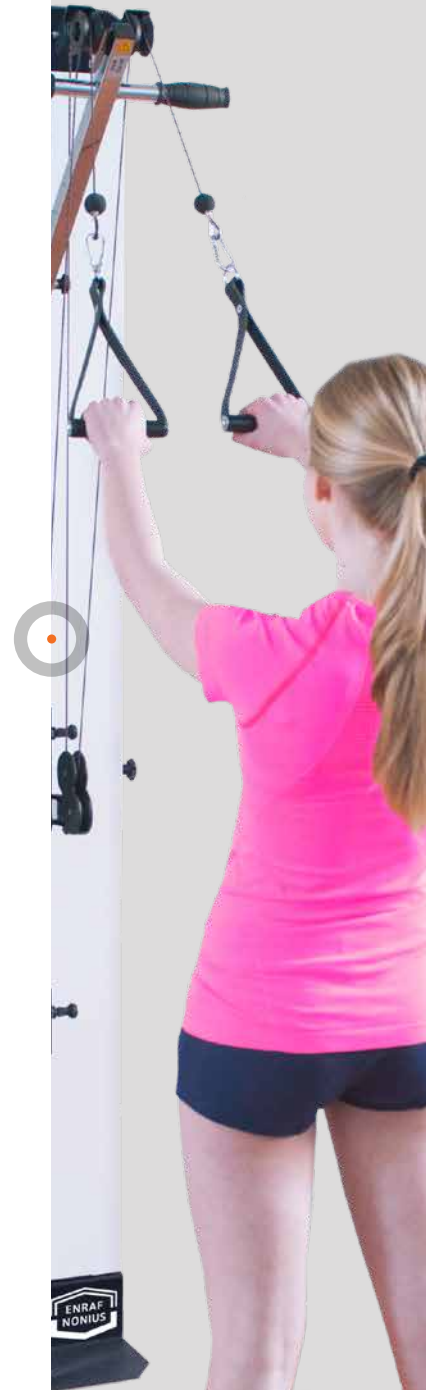

## EN-Tree Pulley Explosive (54 kg)

## ORDERING DATA

The EN-Tree Pulley Explosive is constructed for high speed, explosive movements, which occur for example during sports activities like tennis and soccer. This allows for high velocity movement without high inertia. The EN-Tree Explosive provides the optimal balance between acceleration and resistance for high speed training. The EN-Tree Explosive is provided with 54 kg of weight, which effectively means a maximum weight of 18 kg.

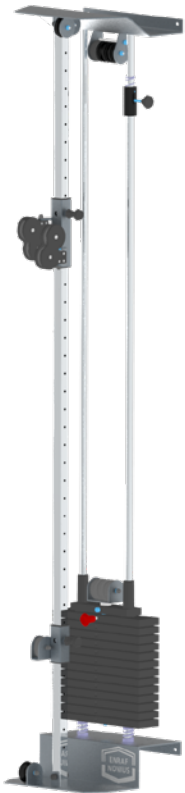

1455972 EN-Tree Pulley Explosive (54 kg)



1455973 EN-Tree Pulley Explosive MD (54 kg, with cover)



# EN-Tree P

# ORDERING DATA

1455974 EN-Tree Pulley (24 kg)



# EN-Tree P

## RELATED PRODUCTS

### EN-Train

The EN-Train concept give you full control over the active rehabilitation process, ensuring optimum treatment at all times. Training schemes can easily be compiled using the training recommendations already preprogrammed in the software.

### Ordering data

1411811 EN-Train installation kit




# EN-Tree M

1455975 EN-Tree M Standard MD (24 kg, with cover)



# ORDERING DATA

1455976 EN-Tree M Explosive MD (54 kg, with cover)

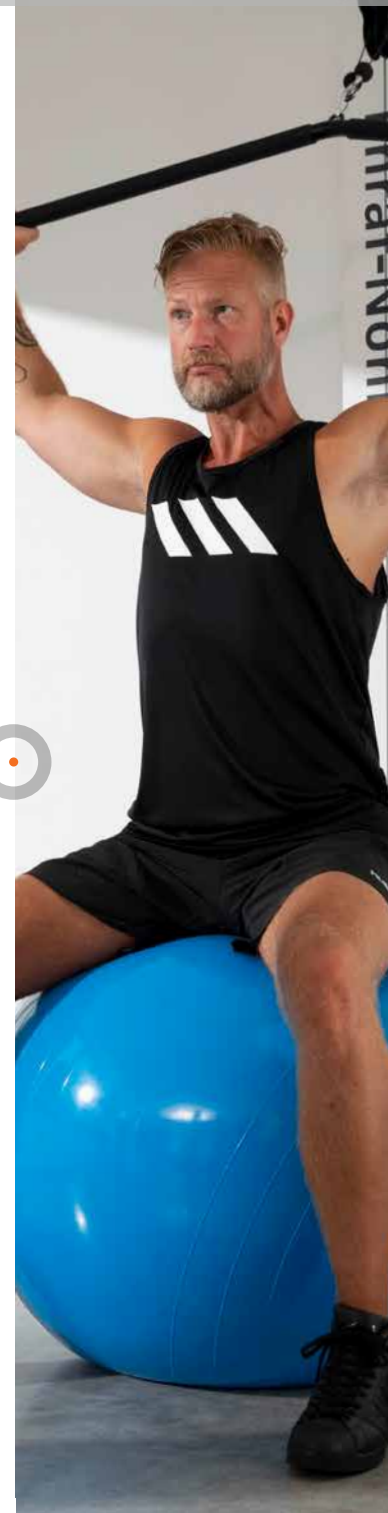

# EN-Tree Bench ORDERING DATA - TECHNICAL SPECIFICATIONS

3443356 EN-Tree Bench

## Two-section MTT training bench

The stable basis for MTT-suitable for both therapeutic support and personal use.

- Practical and easily adjustable
- Inclinable, length-adjustable backrest
- Equipped with inclinometer
- Adjustable seat cushions
- Seat cushions gently inclinable



Dimensions (hwxw):	48x47x134 cm
Weight:	25 kg
Maximum load per section in its lowest position:	135 kg
Angle adjustment:	back-section: 0° to 75° leg-section: 0° and 7°



# EN-Tree

- 3446189 Exercise stool, height 65 cm
- 3443306 Set accessories for EN-Tree Pulley <sup>(1)</sup>
- 3443305 Set of 2 handgrips
- 3443307 Pull down accessory
- 3443308 Lat Pulley handgrip
- 3443309 Triceps handgrip
- 3443310 Ankle strap, 37 cm



## ACCESSORIES



(1) = 3443307 + 3443308 + 3443309 + 3443310

# EN-Tree

# ACCESSORIES

- 3447695 Foot sling
- 3443313 Leather strap, 73 cm
- 3445553 Padded hip belt
- 3445552 Padded ankle belt
- 3443383 Wrist trainer
- 3443381 Bent triceps bar
- 3496253 Row handle



3447695



3443313



3445553



3445552




3443383



3443381



3496253

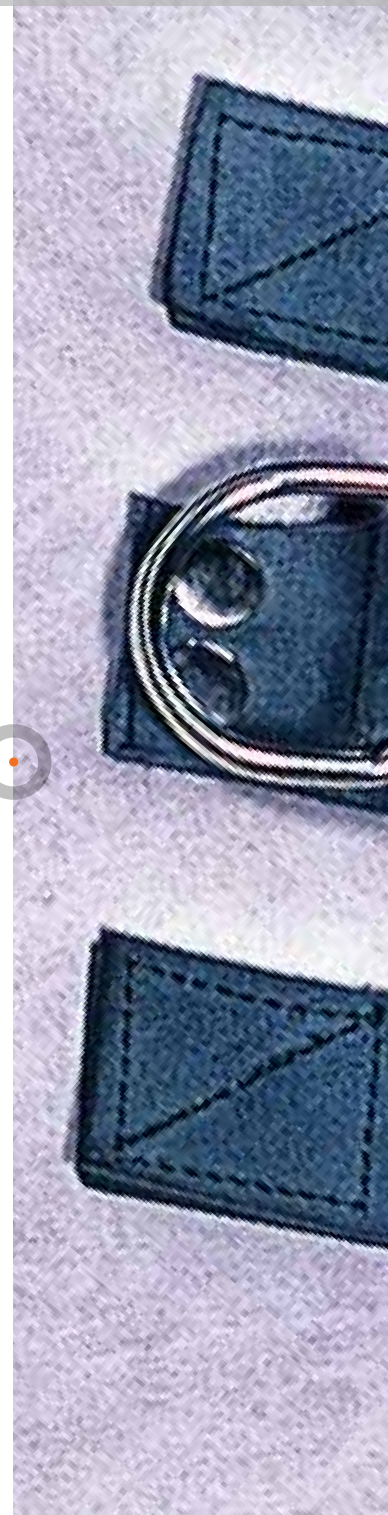

# EN-Tree

# ACCESSORIES

		EN-Tree 24 kg	EN-Tree 54 kg	EN-Tree P	EN-Tree M 24 kg	EN-Tree M 54 kg
3443385	Set of free weights (100 and 200 g)	•			•	
3443386	Set of free weights (300 and 600 g)		•			•
3444008	Mains adaptor 230V/9V - EUR					•
3444032	Mains adaptor 240V/9V - UK					•
3444019	Mains adaptor 115V/9V - US					•



3443385



3443386




3444008



3444032





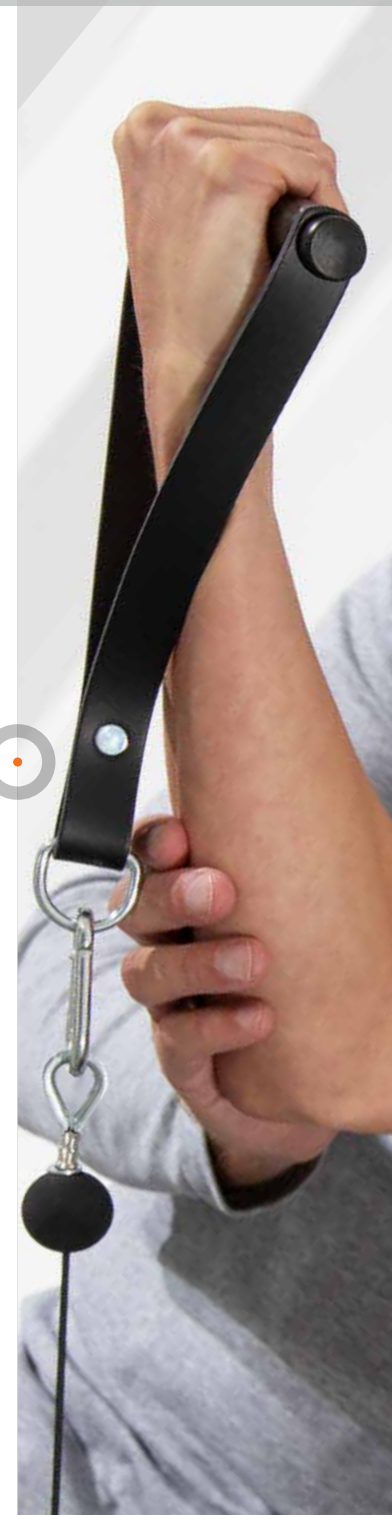




3444019




# EN-Tree

## TECHNICAL SPECIFICATIONS

Articlenumber		Provided with cover	Dimensions (h x w x d) in cm	Weight	Height adjustment for rope	Adjustable weight / load	Effective weight/ load with one cord end	Effective weight/load with 2 cord ends combined	Max. effective weight
1455970	 EN-Tree Pulley (24 kg)		218x33x30	47	0-200	1 kg - 24 kg	A	B	24 kg
1455971	 EN-Tree Pulley MD (24 kg, with cover)	x	218x33x30	60	0-200	1 kg - 24 kg	A	B	24 kg
1455972	 EN-Tree Pulley Explosive (54 kg)		218x39x35	82	0-200	2 kg - 54 kg	C	D	18 kg
1455973	 EN-Tree Pulley Explosive MD (54 kg, with cover)	x	218x39x35	97	0-200	2 kg - 54 kg	C	D	18 kg
1455974	 EN-Tree P	x	218x39x35	60	0-200	40 - 720 N	E	F	24 kg
1455975	 EN-Tree M Standard MD (24 kg, with cover)	x	218x33x30	60	0-200	1 kg - 24 kg	A	B	24 kg
1455976	 EN-Tree M Explosive MD (54 kg, with cover)	x	218x39x35	97	0-200	2 kg - 54 kg	C	D	18 kg

A = 0.25-0.5-0.75-1-2-3-4-5-6-7-8-9-10-11-12 kg  
 B = 0.5-1-1.5-2-4-6-8-10-12-14-16-18-20-22-24 kg  
 C = 0.33-1-1.66-2.33-3-3.66-4.33-5-5.66-6.33-7-7.66-8.33-9 kg

D = 0.66-2-3.33-4.66-6-7.33-8.66-10-11.33-12.66-14-15.33-16.66-18 kg  
 E = 6.66-8.33-10-11.66-13.33 etc ... 120 N (~0.66-0.83-1-1.16-1.33 etc ... 12 kg)  
 F = 13.33-16.66-20-23.33-26.66 etc ... 240 N (~1.33-1.66-2-2.33-2.66 etc ... 24 kg)

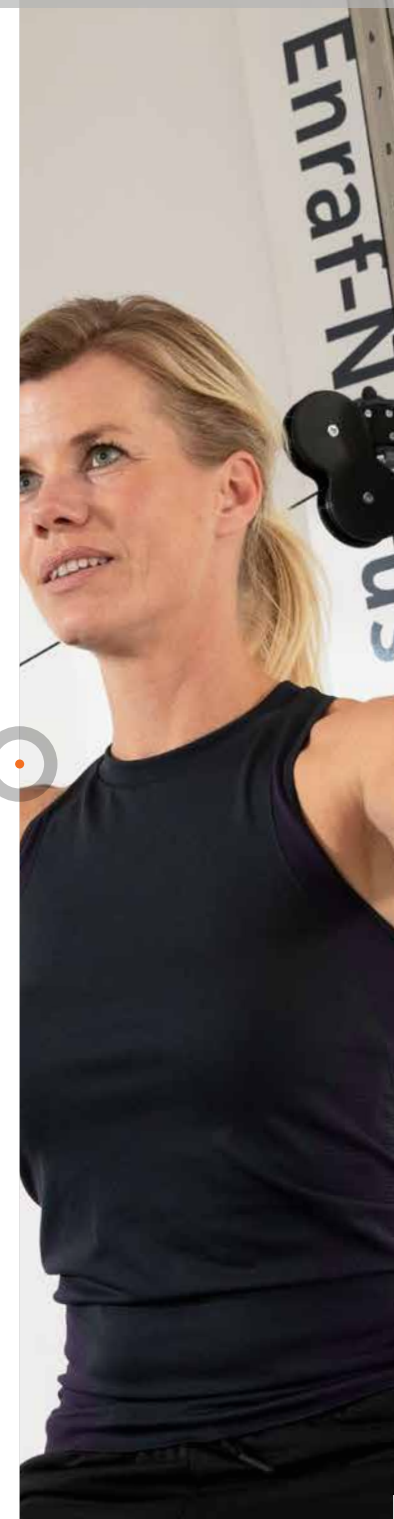

# EN-Tree

## TECHNICAL SPECIFICATIONS

### Installation requirements EN-Tree

EN-Tree should be mounted to the wall with 4 screws and plugs. The wall should be firm enough to resist forces up to 950 N (thus do not attach it to carton-board or thin wooden walls). Leave about 2.5 m<sup>2</sup> free around every pulley for optimal exercise.

### Installation requirements EN-Tree M

The EN-Tree M needs a computer for full operation. The computer can be obtained locally.

### EN-Tree MD

This is a medical device. It complies with all applicable medical regulations. For more details please consult the instructions for use (downloadable via [www.enraf-nonius.com](http://www.enraf-nonius.com)).





## EN-Tree Pulley

A complete series of  
multi-functional pulley systems

Enraf-Nonius B.V. | Vareseweg 127 | 3047 AT Rotterdam | The Netherlands

[www.enraf-nonius.com](http://www.enraf-nonius.com) | [info@enraf-nonius.nl](mailto:info@enraf-nonius.nl) | +31 - (0) 10 20 30 600



PARTNER FOR LIFE